



Barkey S-line

Barkey S-line is a one-touch-control infusion and blood warmer. The device helps to accelerate the healing process of patients and conserves hospital's budgets.

Benefits at a glance

- Warms infusion solutions, FFP and blood up to the patient
- Designed for any standard I.V. set
- One preset temperature for easy operation
- One-touch control
- No additional disposables required
- Displays actual temperature of the warming sleeve

- Use of biocompatible materials
- Visual and acoustic alarm in case of error
- Available in two different diameters
 (5 mm, 7 mm) and lengths (1.2 m, 1.5 m)
- Can be applied for the return of dialyses fluids



P O Box 2541 Southport, QLD 4215 Ph: 1300-786-454 F: 5591 8522 E; sales@questsurgical.com.au ABN 33 118 090 618 ACN 118 090 618




Barkey S-line

Safety

- Integrated temperate sensors
- Overtemperature alarm
- Smooth surfaces, therefore simple to wash and disinfect
- Operational safety thanks to preset target temperature value

Prevention of hypothermia

Warmth is the key to patients' wellbeing. Barkey's inexpensive infusion and blood warmer S-line helps to ease the treatment for all patients and accelerates the healing process. As a result care personnel has less work and treatment costs are reduced.

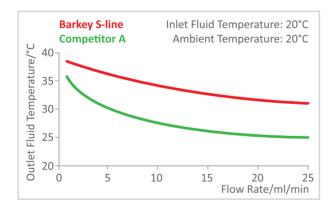



Product	Barkey S-line
Mains supply	100-240 V, 50/60 Hz
Power consumption	115-139 VA
Set temperature	+37°C oder +39°C
Length warming sleeve	1.5 m
For I.V. set diameter	4.0–5.0 mm or 6.0–7.0 mm
Access line	3 m
Weight	approx. 1.3 kg

MDD classification	lla
Protection class	I, protective earthing
Moisture protection	IP X2
Degree of protection	BF, defi-proof /
	internal and external treatment
	to the patient
Certification	CE 0123
Barkey is certified	DIN EN ISO 13485

P O Box 2541 Southport, QLD 4215 Ph: 1300-786-454 F: 5591 8522 E; sales@questsurgical.com.au ABN 33 118 090 618 ACN 118 090 618